

ISO STIM II

Isolated Stimulator

approved electronics - new design





Keep it simple:

Modular approach allows you to:

- choose your desired voltage range and supply source
- choose from various different control options
- choose your output (voltage or current)

**Only pay for what
you really need!**

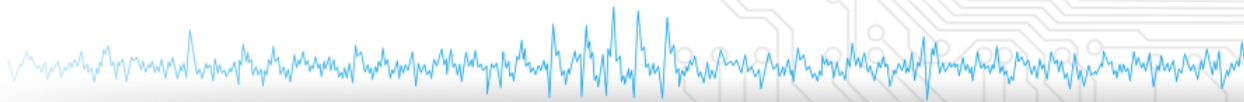

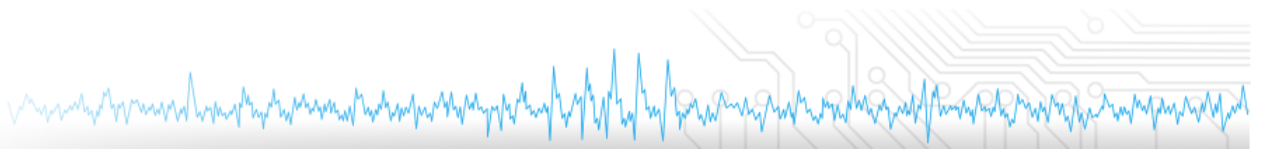

Which options to choose?

Choose base unit with desired voltage range:

ISO-II-100: ± 100 V, ± 10 mA

ISO-II-40: ± 40 V, ± 10 mA

Do you need output other than voltage?

included: voltage output only

ISO-II-CUR: current output instead of voltage output

ISO-II-CURPOT: switchable current **and** voltage output


Choose your control method:

ISO-II-DIR: The output signal follows the input signal




ISO-II-TTL: TTL Gated: Amplitude set at front panel switch

Duration defined by TTL input

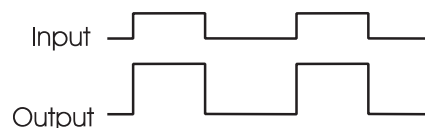

TTL Timed: Amplitude set at front panel switch

Duration set by front panel switch*

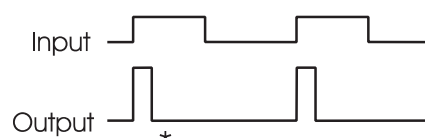




Choose your control method:

All timed pulses can be put out as **unipolar** or **bipolar** pulses



unipolar

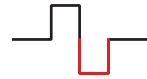


bipolar

*For **bipolar** pulses a variable delay (**IPG**) can be introduced:

ISO-II-IPG: Same features as ISO-II-TTL, **plus:**


InterPhaseGap

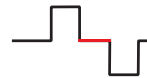

with delay
(interphase gap)

ISO-II-ATR: Duration set by front panel switch*

(Analog Timed) Amplitude follows the input signal

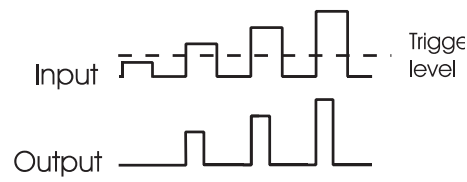

ISO-II-FULL: Includes all above mentioned features

Accessories:

ISO-CAB-D: Shielded cable for electrode connection



DigiStim: Programmable USB signal generator



High quality npi products - made in Germany

Universal Amplifiers

patch clamp, single cell
electroporation, transfection
intra- and extracellular voltammetry

Bridge Amplifiers

most accurate current clamp
amplifiers

Extracellular Amplifiers

LFP, EEG, EMG,
single- or multi-channel

Temperature Controllers

analog control for
low-noise operation

Iontophoretic Substance Application:

fast and precise application of
charged substances

Customized Devices and Prototypes

easy customization due to
inhouse production

Voltammetric Amplifiers

voltammetry, amperometry,
FCV

Miniature Headstages

record in freely moving animals
intra- and extracellular
voltammetry

Amplifiers for Oocytes

fast and accurate clamp with
large currents

Single Electrode Amplifiers

fastest discontinuous current
and voltage clamp

Isolated Stimulators

current and/or voltage
battery or power supply

Filters and Signal Conditioners

high- and lowpass
Bessel filters

Pressurized Substance Application

precise application of
non-polar agents

Fiber Photometry Systems

plug-and-play addition to any
e-phys rig

Contact:

General:
npi electronic GmbH
Phone: +49-7141-9730230
Fax: +49-7141-9730240
sales@npielectronic.com
www.npielectronic.com

North America:
ALA Scientific Instruments
Phone: +1-631-393-6401
Fax: +1-631-393-6407
sales@alascience.com
www.alascience.com